
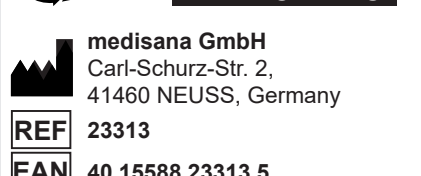
ecomed[®]

by medisana

HM-20E

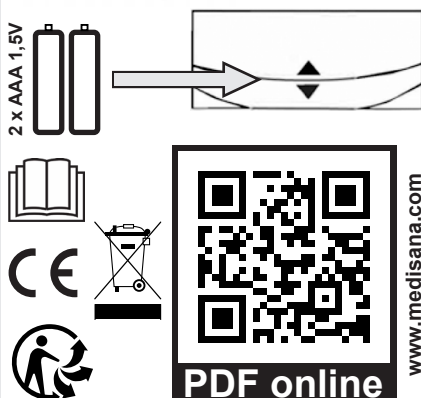www.medisana.com

Mini-Handmassagegerät HM-20E

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Mini Handheld Massager HM-20E

- 1 ON-/OFF-button
 - 2 Hooking nose for massage pad attachment (totally 4 noses)
 - 3 Massage heads (totally 3 massage heads)
- Safety Instructions**
Read the instruction manual carefully before using this device, especially the safety instructions, and keep the instruction manual for future use. These brief instructions are an integral part of the instructions for use. The complete instructions for use are available for download at <https://www.medisana.com>. Should you give this device to another person, it is vital that you also pass on the instructions for use.

- Please ensure that the packing is kept away from the reach of children! Risk of suffocation!
- The device is intended for private use only. If you have any health concerns, consult your doctor before use.
- The device must only be used for its intended purpose as described in the instruction manual. In case of misuse, the warranty claim expires.
- If you experience any pain or discomfort during the massage, stop immediately and consult your doctor.
- Do not use the unit if you suffer from one or more of the following illnesses and complaints: Circulatory disorders, varicose veins, open wounds, bruises, skin tears, phlebitis or thrombosis.
- Do not use the unit if you are pregnant. Special care is necessary when the appliance is used on, by or near children, sick and helpless persons.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.
- Do not use this unit to supplement or replace medical treatment. Chronic complaints and symptoms may be made worse.
- Do not use the appliance near your eyes or near any other sensitive areas.
- Very long periods of use can lead to overheating. So after continuously operating for 15 minutes, user should stop it for at least 15 minutes before using it again.
- Do not use the appliance if it shows signs of damage, if it is not in perfect working order, or if it has been dropped or dipped into water. In order to prevent hazards, always send the unit in to the service centre.
- Only use accessories intended for the appliance and supplied by the manufacturer.
- Do not immerse the device in water or other fluids.
- Do not insert objects into the openings of the device.
- Do not operate the unit where aerosol (spray) products are being used or where oxygen is being administered.
- If the unit is not going to be used for a long period, please remove the batteries.

Battery Safety Information

- Do not disassemble batteries!
- Clean the battery / item contact points before inserting the batteries if necessary!
- Remove discharged batteries from the device immediately!
- If a battery has been swallowed seek medical attention immediately!
- Increased risk of leakage, avoid contact with skin, eyes and mucous membranes! If battery acid comes in contact with any of this parts, rinse the

- affected area with copious amounts of fresh water and seek medical attention immediately!
- Replace all of the batteries simultaneously!
- Only replace with batteries of the same type, never use different types of batteries together or used batteries with new ones!
- Insert the batteries correctly, observing the polarity!
- Remove the batteries from the device if it is not going to be used for an extended period!
- Keep batteries out of children's reach!
- Do not attempt to recharge these batteries! There is a danger of explosion!
- Do not short circuit! Do not throw into a fire! Danger of explosion!
- Do not throw used batteries into the household refuse; put them in a hazardous waste container or take them to a battery collection point, at the shop where they were purchased.

Power supply: 3,0 V⁻⁻⁻, 2 x 1,5V batteries, type AAA (LR03)

Disposal: This product should not be disposed of together with domestic waste. Please hand in all electrical or electronic devices at a municipal or commercial collection point so that they can be disposed of in an environmentally acceptable manner. Consult your municipal authority for information about disposal. Please remove the batteries before disposing of the equipment. Do not dispose of old batteries with your household waste, but dispose of them at a battery collection station at a recycling site or in a shop. Consult your municipal authority or your dealer for information about disposal.

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